

LT5 — Receive DICOM files

1. Open and Configure the DICOM Receiver.


- Go to C:\LT5 folder
- Right click on file “UDTCPDR1.exe”
- Click on Run as administrator.

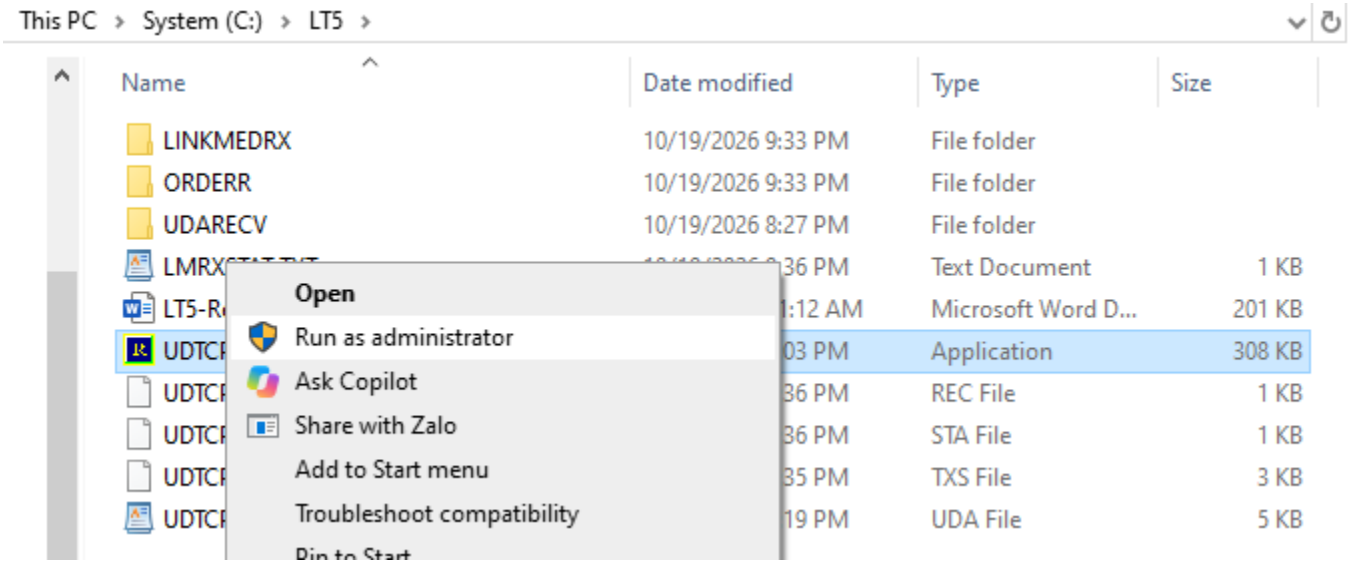


- If Windows displays a security warning, click Yes or Allow, as appropriate.
- The Receiver main window will appear as shown below:

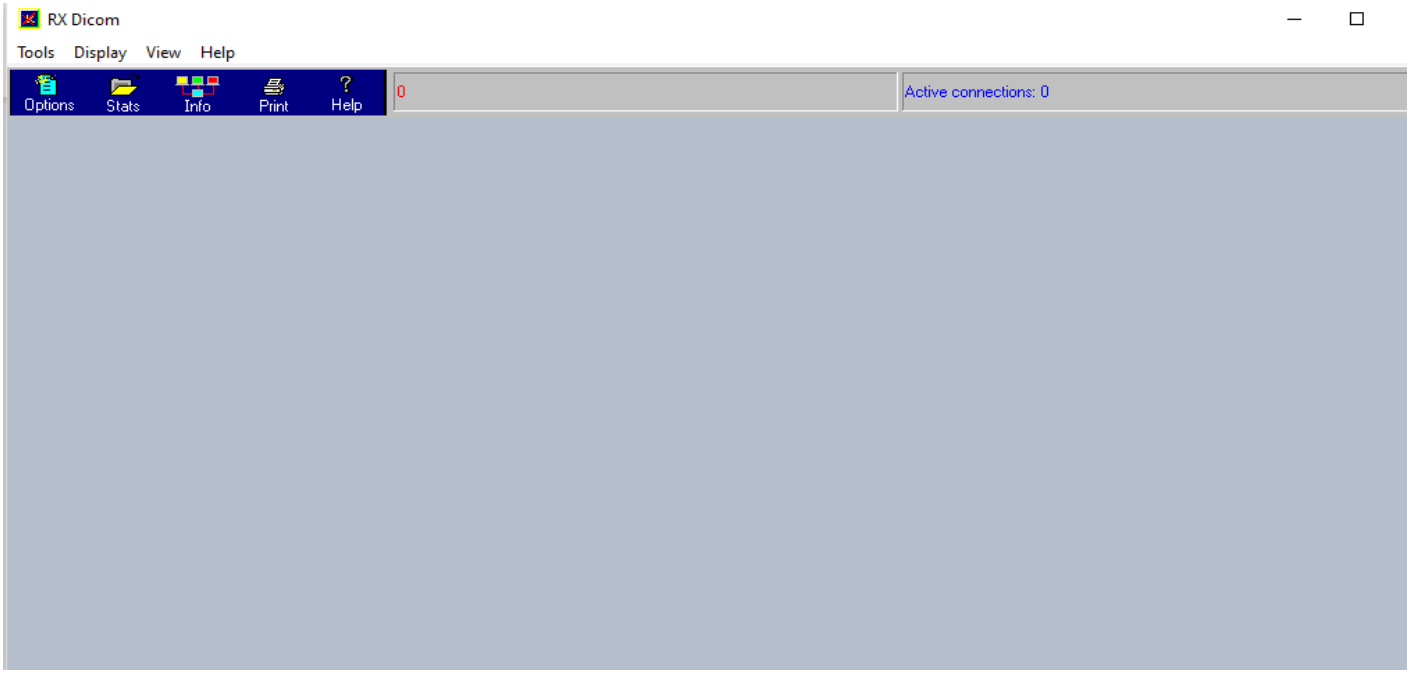


- Next, click on “**Option**” button.

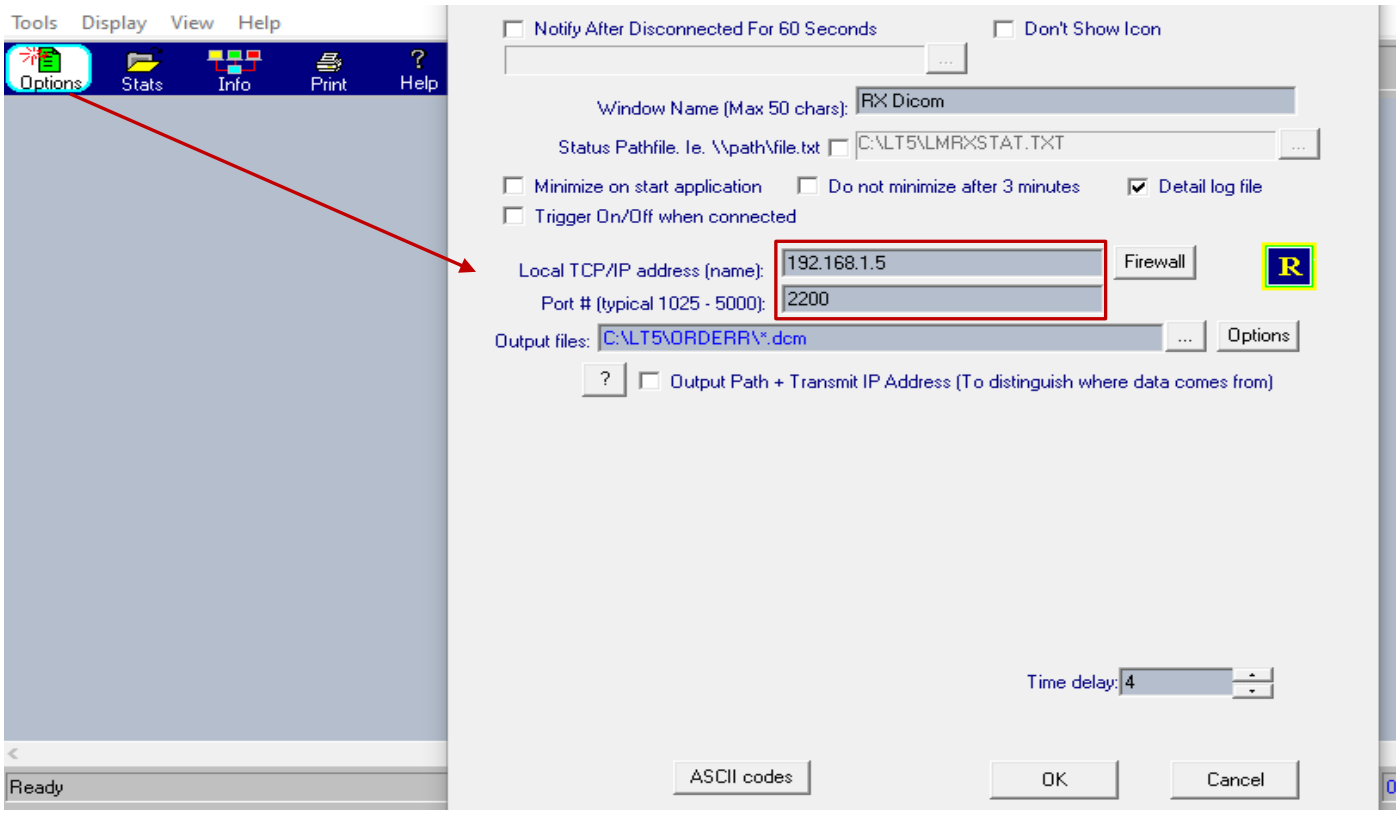


- The Local TCP/IP address and Port # are the IP address and port number of your computer. Provide these to the DICOM sender so it knows where to send the files.
- Change the Local TCP/IP address to the IP address of your computer.
- Leave the default port # 2200.
- Click OK to complete the Receiver configurations.
- When prompted to restart the application, click Yes.

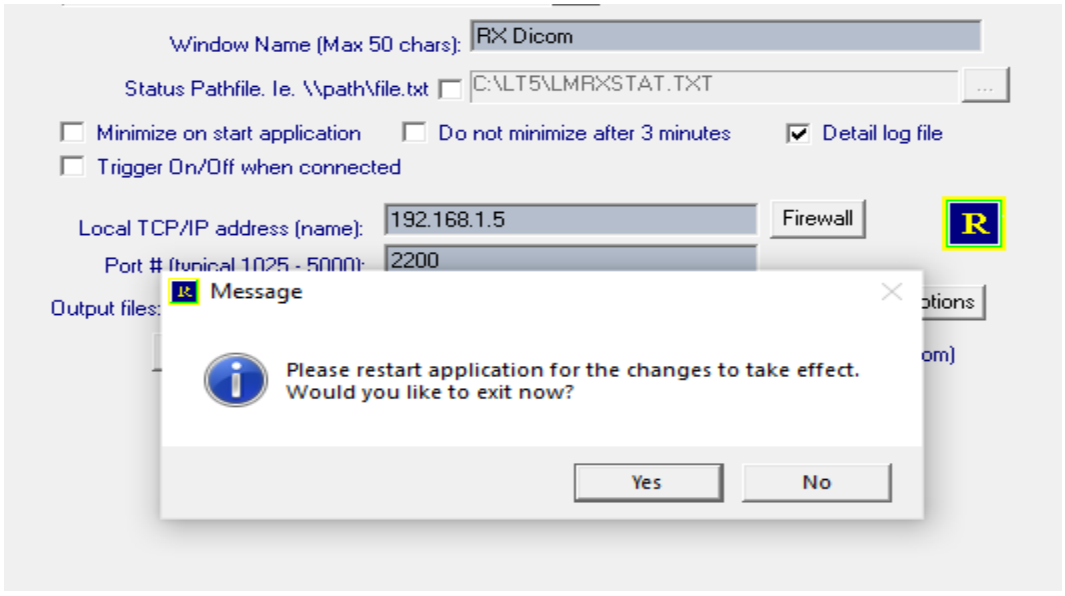

- When asked to confirm that you want to exit, click Yes.

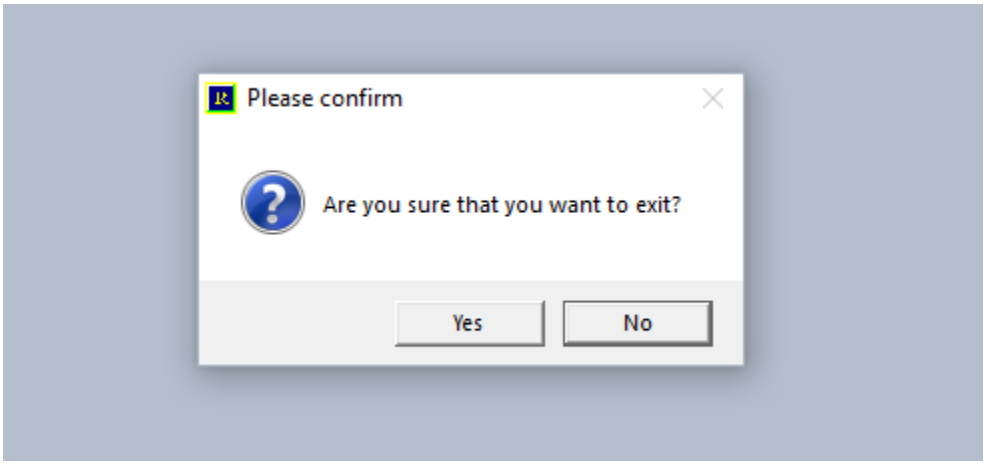

- Go back to the “C:\LT5” folder. Right-click “UDTCPDR1.exe” and select Run as administrator to restart the Receiver with the new configuration.

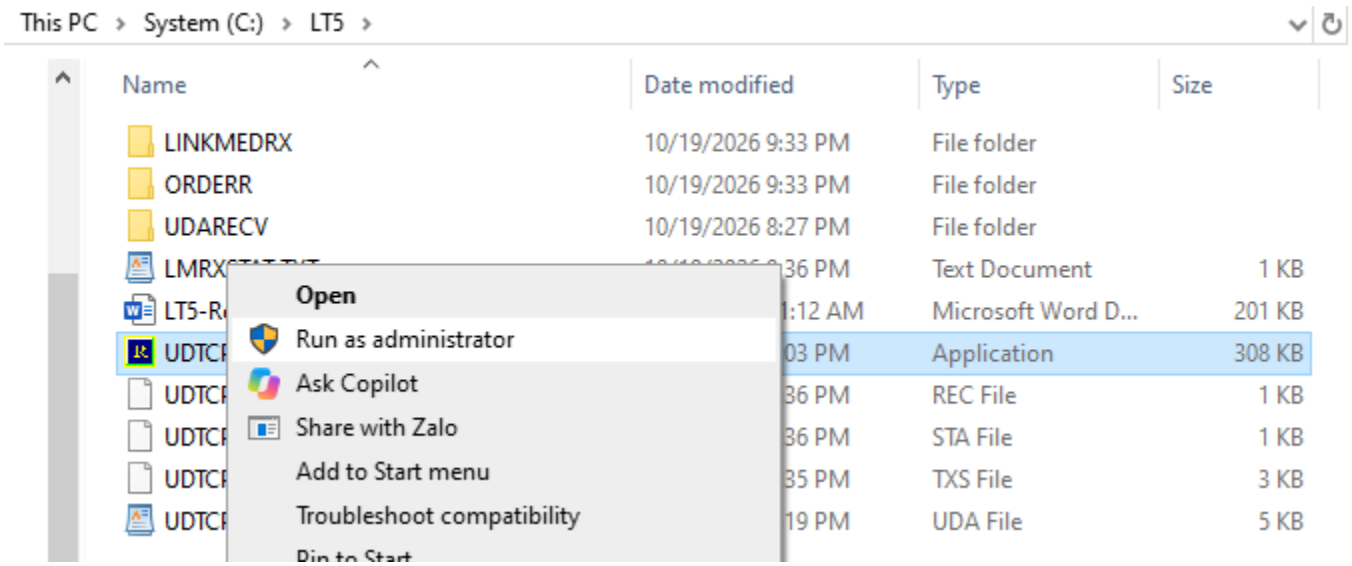

- The Receiver is now listening on your local PC using your local IP address and port 2200.
- If you do not have a DICOM sender for testing, see linkmed.com/linktools/lt6/ to download LT6 and follow the instructions to send a test DICOM file.